
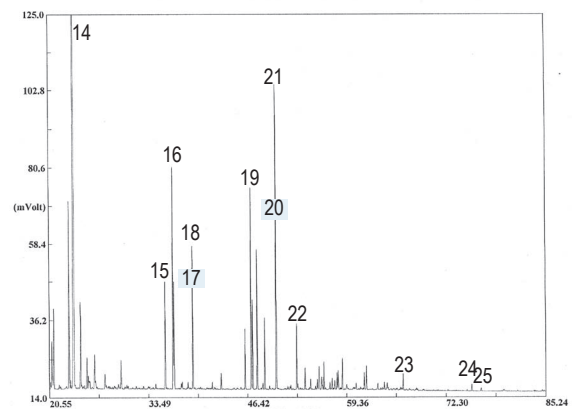


PETROL

Column: **TRB-PETROL**, P/N TR-110592
 Dimensions: 100m x 0.25mm x 0.50µm
 Injection: 0.1µL petrol, split 100:1, 280°C
 Carrier gas: H₂, presión constante 221 Kpa (35°C)
 Oven temperature: 35°C(18min) @ 2°C/min to 200°C(5min)
 Detector: FID, 280°C

Peak Name

1	isobutane	13	n-heptane
2	n-butane	14	toluene
3	isopentane	15	ethylbenzene
4	pentane	16	m-xylene
5	2,3-dimethylbutane	17	p-xylene
6	2-methylpentane	18	o-xylene
7	3-methylpentane	19	1-methyl-3-ethylbenzene
8	hexane	20	1,3,5-trimethylbenzene
9	2,4-dimethylpentane	21	1,2,4-trimethylbenzene
10	benzene	22	1,2,3-trimethylbenzene
11	2-methylhexane	23	naphtalene
12	3-methylhexane	24	2-methylnaphtalene
		25	1-methylnaphtalene

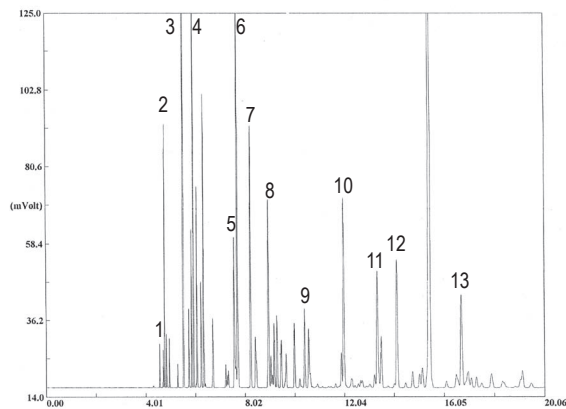

TKG 1203

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es